

SCOPE HELP**O.R.**™



The smart way to achieve clear, surgical vision


The **Scope HelpOR™** is a smart, cost-effective scope management system that helps to optimize surgeons' vision during traditional and robotic assisted laparoscopic surgery.

By combining the tools that reduce condensation and scope fogging with those needed to keep the scope lens and trocar valves free of debris, the **Scope HelpOR™** provides everything needed for effective scope management.

The **Scope HelpOR™** reduces the need to cleanse the scope lens throughout the surgery, which results in improved visibility, less surgical interruption, and therefore fewer procedure delays.



Lap Scope WarmOR™

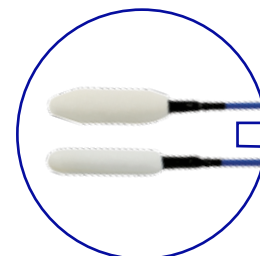

Trocar swabs

CONTACT US

Asia Pacific
122 Balmain Street
Richmond, VIC, 3121
Australia

T +61 3 9413 5555
F +61 3 9413 5556

North America
1625 South Tacoma Way
Tacoma, WA, 98409
USA

T +1 253 441 6534
F +1 866 333 9524

info@theorcompany.com
www.theorcompany.com

Warms and defogs scopes

- Disposable, sterile scope warmer
- Warms the laparoscopes to body temperature in 5-10 minutes, reducing the condensation that causes the scope to fog up during the initial 15-20 minutes of the procedure
- Plug & Play application
- Stays warm throughout the entire procedure
- Accommodates scopes up to 12mm in size, including bariatric and robotic scopes
- Medical grade, reusable 6V AC transformer*



Without Scope HelpOR™ - Fogged view

Cleans scopes and trocars

- Microfiber cloth wipes away persistent smudges without leaving any lint, dust or fibers behind
- Trocar swabs, available in 5-8mm and 10-12mm size, are inserted into the cannula to remove trocar valve debris
- Trocar swab shafts are radiopaque, their rigid, yet flexible, design allows for effective cleaning in tight spaces
- Soft, non-abrasive surgical foam tip absorbs any fluid that has built up



With Scope HelpOR™ - Clear surgical view

* To be purchased separately

Ordering Information

REF/Order #	Description	UOM
SHK-01	Scope HelpOR PRO™ – Lap Scope WarmOR™ warming sheath, trocar swabs (5-8mm and 10-12mm), microfiber cloth	Box of 20, sterile
SW-PA-6V	Lap Scope WarmOR™ – Medical grade, reusable 6V AC transformer	Each, Non-sterile
NLM-3000	Lap Scope WarmOR™ – 7.4V reusable Lithium battery	Each, Non-sterile

CONTACT US

Asia Pacific
122 Balmain Street
Richmond, VIC, 3121
Australia

T +61 3 9413 5555
F +61 3 9413 5556

North America
1625 South Tacoma Way
Tacoma, WA, 98409
USA

T +1 253 441 6534
F +1 866 333 9524

info@theorcompany.com
www.theorcompany.com

ABOUT US

The O.R. Company develops, supplies and markets high quality, innovative surgical devices - from niche products to proprietary surgical instruments - for minimally invasive and open surgery.

Through collaboration with surgeons, clinical technicians and nurses around the world, we identify surgical needs and advance them into effective and safe products such as the Scope HelpOR™.